

THE SYSTEM 2450™





The System 2450[™], a member of ConMed's family of electrosurgical generators, provides state-of-the-art innovation in an easy-to-use format. The System 2450 is simple on the outside, yet sophisticated and powerful on the inside.


Dynamic Response Technology

The System 2450 incorporates ConMed's Dynamic Response technology, the most advanced generator design on the market. Dynamic Response reacts to changes in tissue in less than 10 milliseconds ensuring the delivery of the desired amount of power while limiting the chances of thermal damage to surrounding tissue. Generators without Dynamic Response can produce more arcing and potentially more damage to the surrounding area.

Safe

Patient safety is always a priority. The System 2450 ESU features ConMed's Automatic Return Monitor (ARM) technology designed to alert the operating room staff if the patient/grounding pad contact is compromised when used with compatible dual pad dispersive electrodes. ARM helps to reduce the risk of pad site burns.

When a properly functioning single dispersive electrode is connected to The System 2450, the ARM verifies the connections between the ESU, the dispersive electrode cable, and the dispersive electrode.

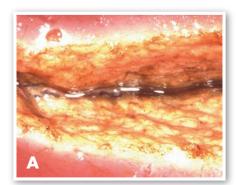

Comparison of Pure Cut at 50 watts on bovine liver. A. ConMed's System 2450 ESU. B. Competitor's ESU. Less charring occurs with The System 2450's Dynamic Response technology.

Smooth Cutting

The System 2450 offers smooth cutting through different tissue types, offering both:

- » Cut for scalpel-like cuts
- >> **Blend** for hemostasis while cutting

Powerful Coagulation

The System 2450 provides powerful coagulation for both dissection and hemostasis. Choose either:

- >> Standard for dissection with hemostasis
- >> **Spray** for broader hemostasis

Precise Bipolar Dessication

The System 2450 provides efficient dessication with minimal lateral thermal tissue damage. Choose:


- » Micro for fine and precise coagulation
- » **Macro** for more energy-intensive coagulation



Features

- » Dynamic Response Technology: The most advanced generator design available
- » Continuous ARM safety monitoring
- » Automatic programming restores The System 2450 to the last settings used
- » One handswitch/footswitch combination jack and a separate handswitch jack allows for multiple accessory connections

- » Bovie-style footswitch connector incorporated into the monopolar handswitch/footswitch jack
- » A serial interface port enables advanced communications between The System 2450 and compatible robotic surgical systems or O.R. integration systems
- >> Integrated interface for activation of a smoke evacuator
- » Front panel volume control

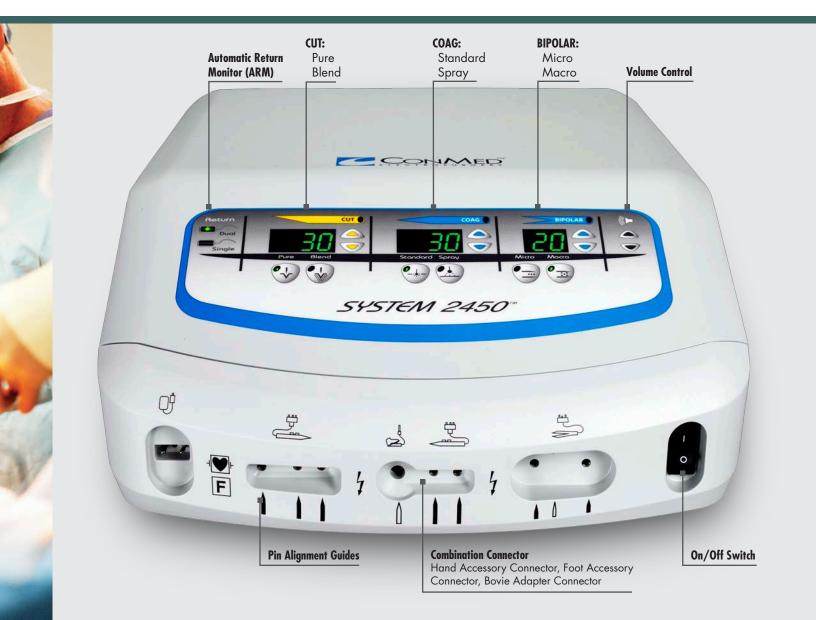

Product Specifications

Watts	Max Peak Voltage	Crest Factor (Typical)	Frequency (KHz)
300	1100	1.4 - 2.0	390 - 410
200	1200	2.0 - 2.6	390 - 410
200	1300	2.2 - 2.9	390 - 410
120	3200	5.6 - 6.6	485 - 505
80	5000	7.5 - 9.5	485 - 505
70	150	1.4 - 2.0	390 - 410
70	450	1.4 - 2.0	390 - 410
	300 200 200 120 80 70	Watts Peak Voltage 300 1100 200 1200 200 1300 120 3200 80 5000 700 150	WattsPeak Voltage(Typical)30011001.4-2.020012002.0-2.620013002.2-2.912032005.6-6.68050007.5-9.5701501.4-2.0

Output Control System	Energy Synchronous Processing
Previous Setting Memory	Yes
Warranty	2 Yrs Standard - 7 Yrs Extended
Contact Quality Monitoring System (CQM)	ARM
CQM Range	10-150 Ohms
Cooling System	Convection
Self Diagnosis	Yes
Technology	Isolated RF Output
Duty Cycle	15 Seconds On, 30 Seconds Off
Compliant with applicable IEC 6060	01 series standards

D

F

11

Weight 32 lbs (14.5 kg)

H: 7.0" (17.8 cm) Dimensions W: 14.4" (36.6 cm) D: 21.5" (54.6 cm)

Please contact your local CONMED Representative for information on our family of generators and our accredited training programs.

Technical Specifications subject to change

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY

1 1

CE

CONMED

e. . SYSTEM 2450



525 French Road Utica, NY 13502 (315) 797-8375 Customer Service: 1-800-448-6506 FAX: 1-800-438-3051 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235 email: info@mail.conmed.com www.conmed.com

© CONMED Corporation, 01/06, printed in U.S.A., Control #12835

Ordering Information

Product No.	Product Description
60-2450-100	System 2450™ ESU (100V)
60-2450-120	System 2450™ ESU (120V)
60-2450-220	System 2450™ ESU (220V)
60-2450-230	System 2450™ ESU (230-240V)
60-7152-001	System 2450 [™] Cart
60-6963-001	Stacking Adapter for System 2450™
60-6700-001	Monopolar Footswitch with 15 ft. cable
60-5103-002	Bipolar Footswitch with 10 ft. cable
60-2451-ENG	Operator's Manual, System 2450 [™] ESU
60-2453-001	In-Service CD, System 2450 [™] ESU
60-2454-ENG	Service Manual, System 2450 [™] ESU
60-8050-001	System 1200 [™] Smoke Evacuator, 120V
60-8055-001	System 1200 [™] SmartFilter [™] VLSI Grade ULPA rated (2 each)

CONMED

(Ų

€₽

JQA

MBR9093

9